


Textile implants for

# BREAST RECONSTRUCTION



Textile  
implantats

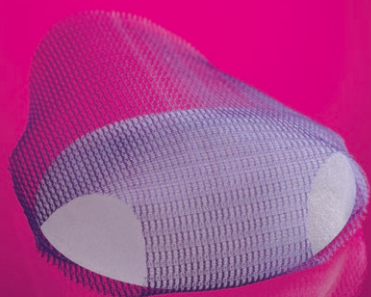


100%  
polydioxanone

completely  
absorbed after  
about 30 weeks

**NEW!** SERASYNTH®  
MESH BR Pocket

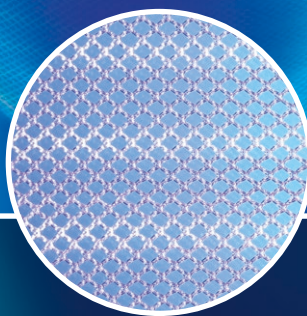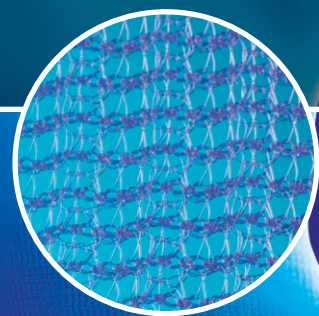


**SERAG**  
**WIESSNER**

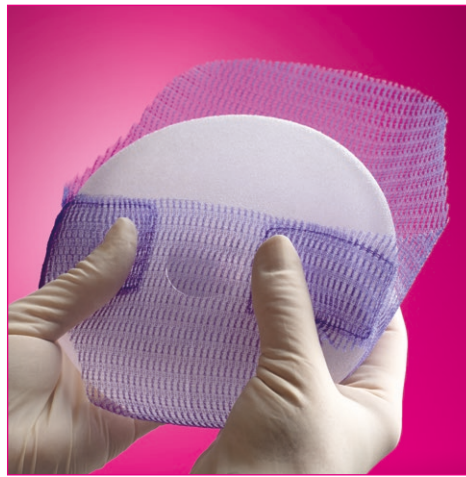

The SERASYNTH® MESH BR Pocket is an absorbable surgical implant for the temporary support of the body's own structures. The ready-made mesh pocket can be used in a pre-pectoral, implant-based breast reconstruction. The pocket can be affixed directly to the pectoralis major muscle or to its fascia so that the muscle

does not have to be detached from the thorax. The front of the implant, which faces the skin, is completely covered with mesh material, while the back is only about half covered.

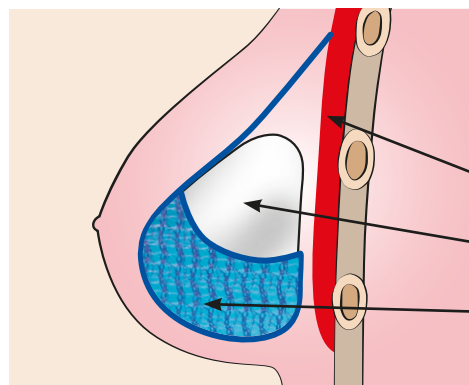

### Pre-pectoral implant placement

Pectoralis major muscle

Silicone implant

SERASYNTH® MESH BR Pocket

Symbol	Name	Packaging Unit	Pore size	Dimension / size*	Order number
	SERASYNTH® MESH BR Pocket	1	~ 1 x 2 mm	14.5 x 17 cm size S ≤ 280 ml	SN44P150
		1	~ 1 x 2 mm	16 x 18.5 cm size M ≤ 420 ml	SN44P250
		1	~ 1 x 2 mm	18 x 20.5 cm size L ≤ 550 ml	SN44P350

\* the sizes of the silicone implants/prostheses (not included in the scope of delivery) are guide values

Risks to be expected in connection with the intended use of the products are specified in the instructions for use. These are available on the SERAG-WIESSNER homepage at [www.serag-wiessner.de](http://www.serag-wiessner.de).



SERAG-WIESSNER GmbH & Co. KG  
Zum Kugelfang 8-12  
95119 Naila/Germany

+ 49 9282 937-230  
+ 49 9282 937-9369  
info@serag-wiessner.de  
www.serag-wiessner.de